



Lightly Powdered Non Sterile Latex Examination Gloves

Kaltex lightly powdered non sterile latex examination gloves with peppermint fragrance are made from Natural Rubber Latex. These mint coloured and flavoured gloves are specially designed for dental use.

FEATURES


- ▶ Non-Sterile, Disposable
- ▶ Specially designed for Dental use
- ▶ The mild peppermint fragrance/flavour is appealing to the patients
- ▶ Unique and special honey comb texture for enhanced grip on instruments
- ▶ Enhanced elasticity, tear strength and donning.
- ▶ Beaded cuff with minimum length of 240 mm.
- ▶ Powdered with Modified Bio-Absorbable Corn starch, to prevent wound infections and granulomas.

STANDARDS

- ▶ EN 455-1: 2020, EN 455-2:2015, EN 455-3:2015 & ASTM D-3578: 2019.

QUALITY CONTROL

- ▶ Factory standard of 1.5 AQL

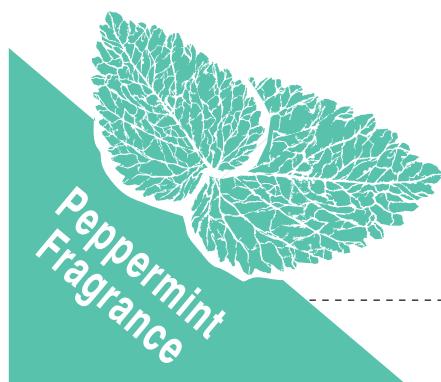

www.kanamlatex.com



Kanam
KANAM LATEX
INDUSTRIES PVT. LTD.

The trusted name the world over

Ooppotttil Buildings, K.K. Road
Kottayam-686 001, Kerala, India
Tel: 0481-2560108, 2563513

Fax: 0481-2563614

E-mail: customercare@kanamlatex.com
ho@kanamlatex.com

Product Information

TYPE

Lightly Powdered Non Sterile Latex Examination Gloves

PRIMARY MATERIAL

Natural Rubber Latex (NR).

COLOUR

Mint Green

DESIGN

Ambidextrous, Straight Finger, Beaded Cuff, Honecomb Textured.

POWDER AMOUNT

<10mg/dm².

SHELF LIFE

3 years in original packaging if stored as per storage instructions.



PHYSICAL DIMENSIONS

Size	Length (mm) Min	Palm Width (mm)	Average weight per pair (gm)	Thickness (mm) Min	
				Palm	Finger
X-Small	240	70 - 79	5.3 - 5.7	0.12	0.15
Small	240	80 - 89	5.7 - 6.1		
Medium	240	90 - 99	6.2 - 6.6		
Large	240	100 - 101	6.7 - 7.1		
X-Large	240	110 - 119	7.3 - 7.7		

PHYSICAL PROPERTIES

Characteristics	Before Ageing	After Ageing 100±2°C for 22±0.3hrs
Tensile Strength (Mpa)	18 min	14 min
Ultimate Elongation (%)	650 min	500 min
Force at break (N) Min.	6 During Shelf Life (N)	6 After Challenging (N)

PRE-SHIPMENT QUALITY INSPECTION

Sampling Procedure: ISO 2859 Part I

Sampling Plan: Single sampling, Normal Inspection

Characteristics	Inspection Level	Method	AQL
Freedom from Holes	G-1	ASTM D 5151:2019/EN 455-1:2020	1.5
Visual attributes – Major Defects	G-1	KL/SP/903	1.5
Visual attributes – Minor Defects	G-1	KL/SP/903	2.5
Physical Dimension	S-2 / N=13	ASTM D 3578: 2019/EN 455-2:2015	4.0
Physical Properties	S-2	ASTM D 3578: 2019	4.0
Force at break	N=13	EN 455-2:2015	N/A
Extractable Protein Content	N=3	ASTM D 5712:2015/EN 455-3: 2015	N/A
Powder Amount	N=2	ASTM D 6124: 2017/ISO 21171: 2006	N/A

PACKAGING

Material	Quantity	Dimension (cm)
Inner Box	100 Pcs.	24(L) x 13(W) x 7(H)
Master Carton	1000 Pcs.	45.5(L) x 34.5(W) x 23(H)



Lightly Powdered Non Sterile Latex Examination Gloves

KLUP/1/22/Rev#03